



DOCOBOAPP™
Peripheral Devices:
Recording temperature

Recording temperature

You will have been provided with a thermometer which looks like one of the thermometers below. Please follow the guidance for the type of thermometer that you have been provided with. Your thermometer may not look exactly the same as the one pictured, but will be very similar and work in the same way.

Thermometers can record temperature in both Degrees Celsius and Degrees Fahrenheit. Please ensure yours is set to record temperature in Degrees Celsius. The display will show °C.

For all thermometers:

- Press the on/off button to activate the thermometer
- The thermometer is ready to use when an indication icon appears on the display. This will vary depending on the brand and type of thermometer.

Tympanic / In Ear Thermometer

- Ensure the ear is free from obstruction and excess ear wax.
- Remove the protective cap and place a new, clean hygiene cap on the thermometer. Do this before every use.
- Place the probe in the ear canal directed towards the opposite ear and press the start button.
- The temperature measurement will appear on the display.
- If the thermometer makes a series of short beeps the hygiene cap may not be fitted correctly. Replace the cap and repeat the process again.
- Write the measurement down and dispose of the used hygiene cap.



Infrared / Forehead Thermometer

- Point the infrared end directly towards your forehead, between the brows and approximately 5cm / 2 Inches away.
- Press the start button which may produce a beep sound. The measurement will be displayed on the screen.
- Write the measurement down as the thermometer will power off automatically after around 10 seconds.



Digital Thermometer

- Place the silver tip under the tongue or in the crease of the arm pit.
- When the °C flashes, the measurement is in progress. When the °C stops flashing, the measurement is complete.
- Write the measurement down as the thermometer will power off automatically.




Recording Temperature on DOCOBOAPP™

- You can enter the temperature measurement by either:

Tapping the grey arrows left and right

Or dragging the grey slider side to side

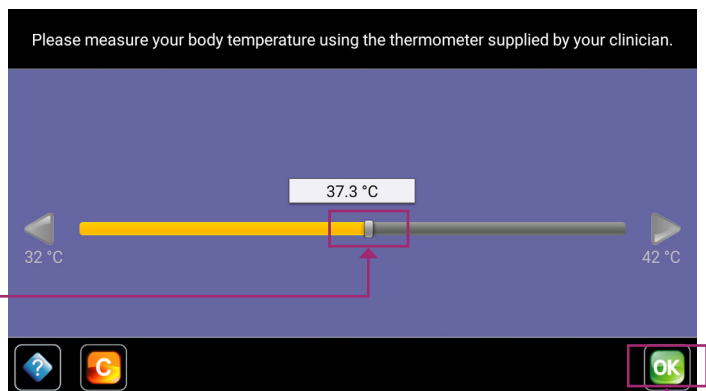

When the measurement is correct press OK



For technical assistance call 01372 899031

Graphnet 
Transforming Care

Docobo is a wholly owned
subsidiary of Graphnet
Health Limited
Docobo
21 High Street Bookham
Surrey KT23 4AA

Graphnet Health Limited
Ground Floor
Building 5 Caldecotte Lake Drive
Caldecotte
Milton Keynes
Buckinghamshire MK7 8LE

T +44 (0)3330 771 988
E info@graphnethealth.com
www.graphnethealth.com